



GE Medical Systems

Technical Publications

Direction 2173227–100

Revision 5

AMX–4+ Periodic Maintenance (Model 2169360, 2236420 & 2275938 Series)

Copyright© 1996, 1997, 1999, 2000 By General Electric Co.

Operating Documentation

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

AVERTISSEMENT

- CE MANUEL DE MAINTENANCE N'EST DISPONIBLE QU'EN ANGLAIS.
- SI LE TECHNICIEN DU CLIENT A BESOIN DE CE MANUEL DANS UNE AUTRE LANGUE QUE L'ANGLAIS, C'EST AU CLIENT QU'IL INCOMBE DE LE FAIRE TRADUIRE.
- NE PAS TENTER D'INTERVENTION SUR LES ÉQUIPEMENTS TANT QUE LE MANUEL SERVICE N'A PAS ÉTÉ CONSULTÉ ET COMPRIS.
- LE NON-RESPECT DE CET AVERTISSEMENT PEUT ENTRAÎNER CHEZ LE TECHNICIEN, L'OPÉRATEUR OU LE PATIENT DES BLESSURES DUES À DES DANGERS ÉLECTRIQUES, MÉCANIQUES OU AUTRES.

WARNUNG

- DIESES KUNDENDIENST-HANDBUCH EXISTIERT NUR IN ENGLISCHER SPRACHE.
- FALLS EIN FREMDER KUNDENDIENST EINE ANDERE SPRACHE BENÖTIGT, IST ES AUFGABE DES KUNDEN FÜR EINE ENTSPRECHENDE ÜBERSETZUNG ZU SORGEN.
- VERSUCHEN SIE NICHT, DAS GERÄT ZU REPARIEREN, BEVOR DIESES KUNDENDIENST-HANDBUCH NICHT ZU RATE GEZOGEN UND VERSTANDEN WURDE.
- WIRD DIESE WARNUNG NICHT BEACHTET, SO KANN ES ZU VERLETZUNGEN DES KUNDENDIENSTTECHNIKERS, DES BEDIENERS ODER DES PATIENTEN DURCH ELEKTRISCHE SCHLÄGE, MECHANISCHE ODER SONSTIGE GEFAHREN KOMMEN.

AVISO

- ESTE MANUAL DE SERVICIO SÓLO EXISTE EN INGLÉS.
- SI ALGÚN PROVEEDOR DE SERVICIOS AJENO A GEMS SOLICITA UN IDIOMA QUE NO SEA EL INGLÉS, ES RESPONSABILIDAD DEL CLIENTE OFRECER UN SERVICIO DE TRADUCCIÓN.
- NO SE DEBERÁ DAR SERVICIO TÉCNICO AL EQUIPO, SIN HABER CONSULTADO Y COMPRENDIDO ESTE MANUAL DE SERVICIO.
- LA NO OBSERVANCIA DEL PRESENTE AVISO PUEDE DAR LUGAR A QUE EL PROVEEDOR DE SERVICIOS, EL OPERADOR O EL PACIENTE SUFRAN LESIONES PROVOCADAS POR CAUSAS ELÉCTRICAS, MECÁNICAS O DE OTRA NATURALEZA.

ATENÇÃO

- ESTE MANUAL DE ASSISTÊNCIA TÉCNICA SÓ SE ENCONTRA DISPONÍVEL EM INGLÊS.
- SE QUALQUER OUTRO SERVIÇO DE ASSISTÊNCIA TÉCNICA, QUE NÃO A GEMS, SOLICITAR ESTES MANUAIS NOUTRO IDIOMA, É DA RESPONSABILIDADE DO CLIENTE FORNECER OS SERVIÇOS DE TRADUÇÃO.
- NÃO TENHA TENTADO REPARAR O EQUIPAMENTO SEM TER CONSULTADO E COMPREENDIDO ESTE MANUAL DE ASSISTÊNCIA TÉCNICA.
- O NÃO CUMPRIMENTO DESTA AVISO PODE POR EM PERIGO A SEGURANÇA DO TÉCNICO, OPERADOR OU PACIENTE DEVIDO A CHOQUES ELÉTRICOS, MECÂNICOS OU OUTROS.

AVVERTENZA

- IL PRESENTE MANUALE DI MANUTENZIONE È DISPONIBILE SOLTANTO IN INGLESE.
- SE UN ADDETTO ALLA MANUTENZIONE ESTERNO ALLA GEMS RICHIEDE IL MANUALE IN UNA LINGUA DIVERSA, IL CLIENTE È TENUTO A PROVVEDERE DIRETTAMENTE ALLA TRADUZIONE.
- SI PROCEDA ALLA MANUTENZIONE DELL'APPARECCHIATURA SOLO DOPO AVER CONSULTATO IL PRESENTE MANUALE ED AVERNE COMPRESO IL CONTENUTO.
- NON TENERE CONTO DELLA PRESENTE AVVERTENZA POTREBBE FAR COMPIERE OPERAZIONI DA CUI DERIVINO LESIONI ALL'ADDETTO ALLA MANUTENZIONE, ALL'UTILIZZATORE ED AL PAZIENTE PER FOLGORAZIONE ELETTRICA, PER URTI MECCANICI OD ALTRI RISCHI.

警告

- ・ このサービスマニュアルには英語版しかありません。
- ・ GEMS以外でサービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
- ・ このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないで下さい。
- ・ この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

注意:

- 本维修手册仅存有英文本。
- 非 GEMS 公司的维修员要求非英文本的维修手册时，客户需自行负责翻译。
- 未详细阅读和完全了解本手册之前，不得进行维修。
- 忽略本注意事项会对维修员，操作员或病人造成触电，机械伤害或其他伤害。

Direction 2173227-100**Revision 5****AMX-4+ Periodic Maintenance
(Model 2169360, 2236420 & 2275938 Series)****IMPORTANT! . . . X-RAY PROTECTION**

X-ray equipment if not properly used may cause injury. Accordingly, the instructions herein contained should be thoroughly read and understood by everyone who will use the equipment before you attempt to place this equipment in operation. The General Electric Company, Medical Systems Group, will be glad to assist and cooperate in placing this equipment in use.

Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical

design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that everyone having anything to do with x-radiation be properly trained and fully acquainted with the recommendations of the National Council on Radiation Protection and Measurements as published in NCRP Reports available from NCRP Publications, 7910 Woodmont Avenue, Room 1016, Bethesda, Maryland 20814, and of the International Commission on Radiation Protec-

tion, and take adequate steps to protect against injury.

The equipment is sold with the understanding that the General Electric Company, Medical Systems Group, its agents, and representatives have no responsibility for injury or damage which may result from improper use of the equipment.

Various protective material and devices are available. It is urged that such materials or devices be used.

CAUTION: United States Federal law restricts this device to use by or on the order of a physician.

THIS PAGE INTENTIONALLY LEFT BLANK.

If you have any comments, suggestions or corrections to the information in this document, please write them down, include the document title and document number, and send them to:

**GENERAL ELECTRIC COMPANY
MEDICAL SYSTEMS**

MANAGER – INFORMATION INTEGRATION, AMERICAS W-622
P.O. BOX 414
MILWAUKEE, WI 53201-0414

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT



All electrical installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations, and testing shall be

performed by qualified GE Medical personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the

requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "**damage in shipment**" written on **all** copies of the freight or express bill **before** delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately**

upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

Call Traffic and Transportation, Milwaukee, WI (414) 827-3449 /

8*285-3449 **immediately** after damage is found. At this time be ready to supply name of carrier, delivery date, consignee name, freight or express bill number, item damaged and extent of damage.

Complete instructions regarding claim procedure are found in Section "S" of the Policy & Procedure Bulletins.

THIS PAGE INTENTIONALLY LEFT BLANK.

TABLE OF CONTENTS

SECTION	TITLE	PAGE
	REVISION HISTORY	ix
1	INTRODUCTION	1-1
	1-1 Identification	1-1
	1-2 General	1-2
	1-3 Cleaning	1-2
	1-4 Inspection	1-3
	1-5 HHS Testing	1-3
	1-6 Tools and Materials	1-3
2	INSPECTION	2-1
APPENDIX		
1	SYMBOLS	A-1

THIS PAGE INTENTIONALLY LEFT BLANK.

REVISION HISTORY

REV	DATE	REASON FOR CHANGE
A	Nov. 6, 1996	Initial draft.
0	Dec. 13, 1996	Initial release.
1	Mar. 7, 1997	Removed reference to 2173467-100.
2	June 24, 1997	Changed inspection frequency to 12 months.
3	Aug. 14, 1997	High Impact Inspection.
4	Apr. 12, 1999	Added AMX-4+ Model 2236420.
5	Nov. 8, 2000	Added AMX-4+ Model 2275938.

LIST OF EFFECTIVE PAGES

PAGE NUMBER	REVISION NUMBER	PAGE NUMBER	REVISION NUMBER	PAGE NUMBER	REVISION NUMBER
Title Page	5				
i thru x	5				
1-1 thru 1-4	5				
2-1 thru 2-2	5				
A-1 thru A-2	5				
Back Page	—				

THIS PAGE INTENTIONALLY LEFT BLANK.

SECTION 1 INTRODUCTION


1-1 Identification

See Illustration 1-1. The AMX-4+ is identified by Model Number on the rating plate located on the top cover. Model part and catalog numbers are identified in Table 1-1.

TABLE 1-1
AMX-4+ MODELS

DESCRIPTION	PART NUMBER	CATALOG NUMBER	PART NUMBER	CATALOG NUMBER
DOMESTIC	2169360-7	A0659F	2236420-7 & 2275938-7	A0659JF
DOMESTIC, AEC	2169360-8	A0659FA	2236420-8 & 2275938-8	A0659JG
DOMESTIC, TECH SWITCH	2169360-9	A0659FC	2236420-9 & 2275938-9	A0659JH
DOMESTIC, AEC, TECH SWITCH	2169360-10	A0659FB	2236420-10 & 2275938-10	A0659JJ
IEC, EMC	2169360	A0659A	2236420 & 2275938	A0659J
IEC, EMC, AEC	2169360-2	A0659AA	2236420-2 & 2275938-2	A0659JA
IEC, EMC, TECH SWITCH	2169360-3	A0659AB	2236420-3 & 2275938-3	A0659JB
IEC, EMC, AEC, TECH SWITCH	2169360-4	A0659AC	2236420-4 & 2275938-4	A0659JC
JAPAN	2169360-5	A0659C	2236420-5 & 2275938-5	A0659JD
JAPAN SHORT COLUMN	2169360-6	A0659D	2236420-6 & 2275938-6	A0659JE

ILLUSTRATION 1-1
AMX-4+ IDENTIFICATION

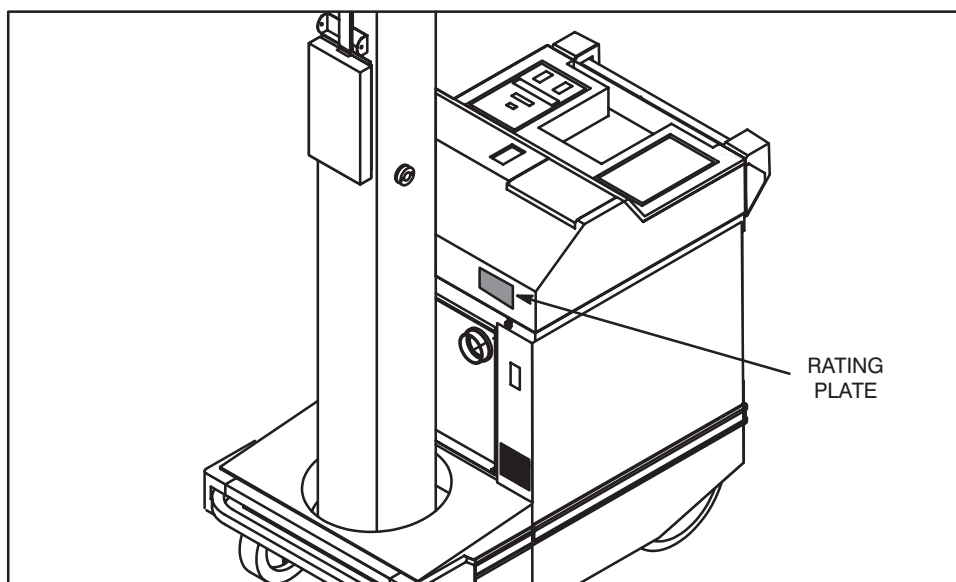

1-2 General

The AMX-4+ contains operating safeguards to provide maximum safety. Before calling for service, be certain proper operating procedures are being used. Refer to Direction 2166913-100, *AMX-4+ Operation*, or Direction 2166911-100, *AMX-4+ International Operation*, for proper operating procedures.

Satisfactory equipment performance requires the use of service personnel specially trained on x-ray apparatus. GE Medical Systems is responsible for the effects on safety, reliability, and performance only if the following conditions are met:

- The electrical wiring of the relevant rooms complies with all national and local codes.
- All assembly operations, extensions, re-adjustments, modifications, or repairs are carried out by GE Medical Systems' authorized service representatives.
- The equipment is used in accordance with the instructions for use. Refer to Direction 2166913-100, *AMX-4+ Operating Manual*, or Direction 2166911-100, *AMX-4+ International Operation*, for proper operating procedures.



Only trained and qualified personnel should be permitted access to the internal parts of this equipment.

GE Medical Systems, and its associates, maintain a worldwide organization of stations from which one may obtain skilled x-ray service. If desired, arrangements can usually be made to furnish preventive and/or emergency service on a contract basis. A GE Medical Systems representative will be glad to discuss this plan.

1-3 Cleaning

The AMX-4+ should be cleaned frequently, particularly if corroding chemicals are present. Use a cloth slightly moistened in warm, soapy water to clean all surfaces. Wipe with a cloth slightly moistened in clean water. Do not use cleaners or solvents of any kind. They may remove markings and could damage the finish and plastic covers.

1-4 Inspection

In order to assure continued performance of this x-ray equipment, a periodic inspection program must be established. It is the owner's responsibility to supply or arrange for this service.

Periodic Maintenance can locate potential problems before they occur. Since maintenance checks involve cost, it is not practical to attempt to assure zero failures. The following inspection requirements, with suggested inspection intervals, are GE Medical Systems' recommendation for the most cost effective maintenance schedule.

Inspection intervals are based on average daily use of one eight hour shift. More frequent inspection is appropriate where equipment use is above average. Equipment that regularly exhibits out of tolerance conditions should also be inspected more frequently.

1-5 HHS Testing

The United States Department of Health and Human Services (HHS) has established performance requirements for diagnostic x-ray equipment. These requirements are defined in Title 21 of the Code of federal regulations (21 CFR), and apply only to certain specified components identified as “certified equipment.”

The manufacturer of specified diagnostic x-ray components must certify that the components:

1. Perform as required by the HHS standard when installed, adjusted, and tested as specified in the manufacturer’s instructions to the assembler.

Refer to Direction 2173222-100 *AMX-4+ Installation* for HHS compliance testing as it pertains to the AMX-4+ during installation.

2. Will continue to comply when maintained in accordance with the manufacturer’s instructions.

The importance of HHS compliance testing cannot be over emphasized. Before starting Periodic Maintenance procedures, review Direction 46-013894 *System Field Test for HHS*. Pay particular attention to “Introduction” (Tab 1 of Direction 46-013894).

1-6 Tools and Materials


In addition to the standard service representatives tool kit, the following items are required:

- Direction 2173221-100, *AMX-4+ Ratings & Specifications* and associated tools and materials.
- Direction 2173224-100, *AMX-4+ Functional Check* and associated tools and materials.
- Direction 2173225-100, *AMX-4+ Service* and associated tools and materials.
- Direction 2190775-100, *Stranded Steel Cable Inspection*, and associated tools and materials

THIS PAGE INTENTIONALLY LEFT BLANK.

SECTION 2
INSPECTION

The recommended interval between maintenance checks is 12 months. The following flowchart shows the recommended sequence for performing an annual periodic maintenance check on the AMX-4+.

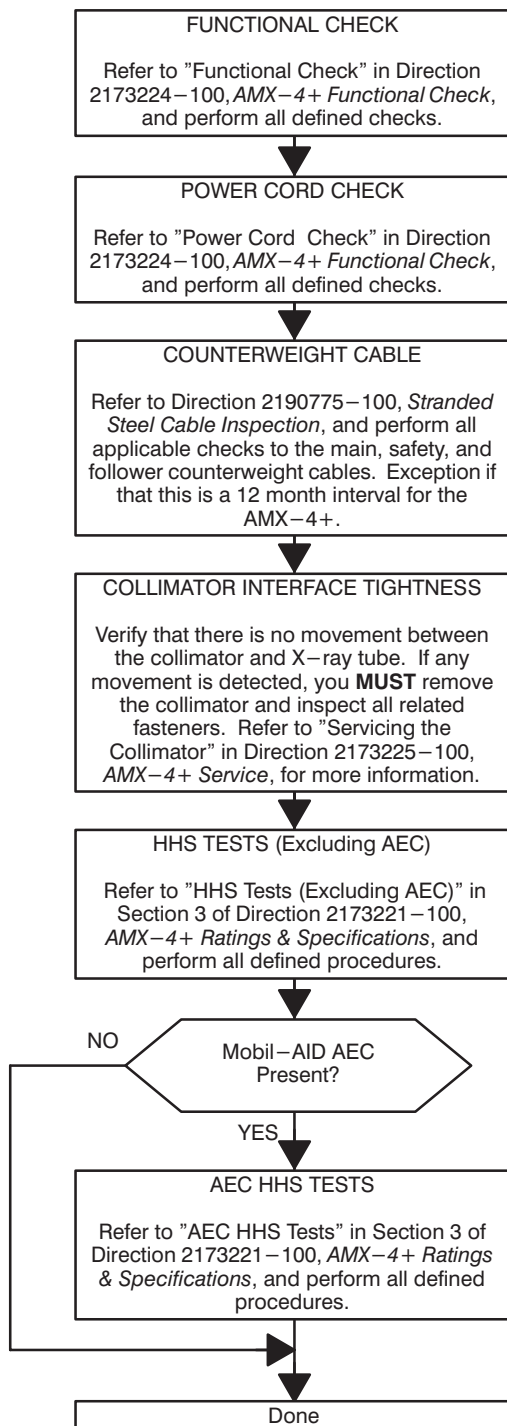ILLUSTRATION 2-1
INSPECTION STEERING

THIS PAGE INTENTIONALLY LEFT BLANK.

APPENDIX 1

SYMBOLS

All symbols used on the equipment and in its accompanying documents are shown and explained in this appendix.



Caution advises of an avoidable condition that could cause minor physical injury, or damage to equipment or data.



Warning advises of an avoidable condition that may allow or cause a personal injury or the catastrophic destruction of equipment or data.



Danger advises of an avoidable condition that will cause serious or fatal injury.



Type B Equipment. Internal electrical power source provides an adequate degree of protection against electrical shock.



X-ray emission. X-ray tube head is emitting x-rays. Take adequate precautions to prevent the possibility of any persons carelessly, unwisely, or unknowingly exposing themselves or others to radiation.



Battery power on. This does not apply mains voltage.



Battery power off. This does not remove mains voltage.



Control for indicating radiation field by using light.



Collimator blades closed. The controlled blades are shown in thicker lines.



Collimator Blades open. The controlled blades are shown in thicker lines.



Functional Earth (ground) Terminal. Terminal directly connected to a point of a measuring supply or control circuit or to a screening part which is intended to be earthed for functional purposes.



Alternating Current. Indicates equipment that is suitable for alternating current only.



Direct Current. Indicates equipment that is suitable for direct current only.



Equipotentiality. Identifies terminals that bring the various parts of equipment or systems to the same potential when connected together. These terminals are not necessarily at earth (ground) potential. The value of the potential may be indicated next to the symbol.



Indicates lock release or brake release.



Indicates receptacle location for hand-held radiographic prep/expose and field-light control cable.



GE Medical Systems

*GE Medical Systems: Telex 3797371
P.O. Box 414, Milwaukee, Wisconsin 53201 U.S.A.
(Asia, Pacific, Latin America, North America)*

*GE Medical Systems — Europe: Telex 698626
283, rue de la Minière, B.P. 34, 78533 Buc Cedex
France*